



TRIFACIAL SERVICE CENTER

Save space, time and money without sacrificing versatility. The **Trifacial Service Center** provides an ample worksurface and central sink allowing the students to easily access the services. When the Trifacial is combined with Sheldon Student Tables, the classroom becomes a lecture and lab environment which reduces classroom disruption and provides more circulation area for the students.

Features


- Three sides of triangular configuration
- Large, 34-gallon capacity sink for easy cleanup
- Center height adjustable at time of installation for elementary, middle school, and high school applications
- Equipped with three vandal-resistant unicast cold water and gas fixtures and three GFI outlets
- Support structure of fire-retardant fiberglass reinforced polyester

sales@sheldonlabs.com



800.531.7604

TRIFACIAL SINK CENTER

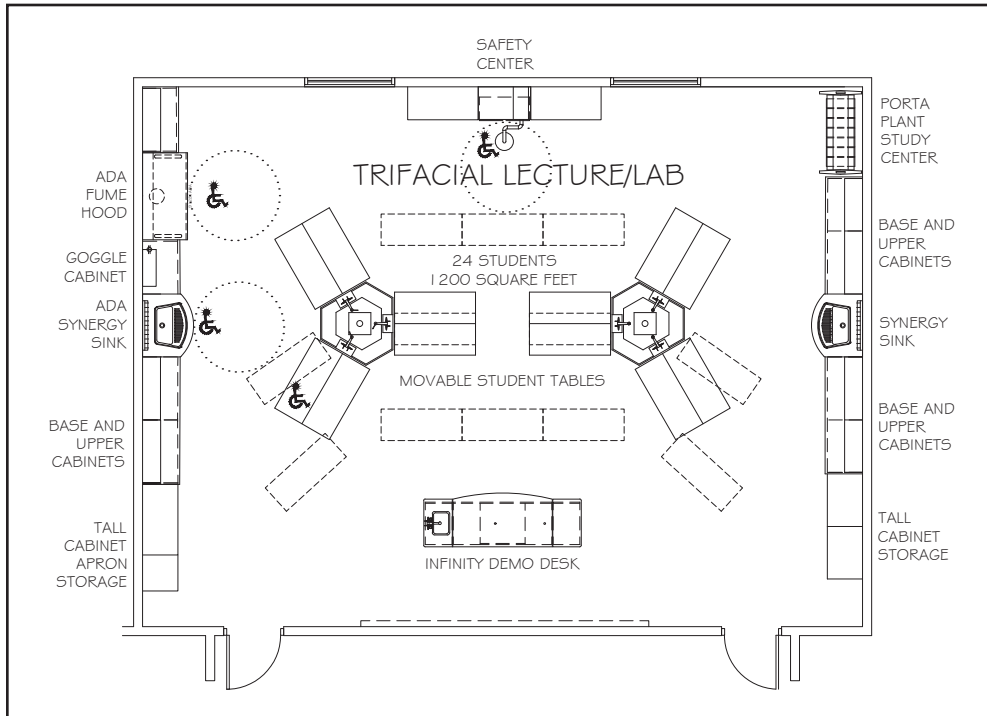

Trifacial Lab

This lab is designed as lecture and lab layout. This allows the students to collaborate with each other while creating a sightline for the instructor to all locations. The student tables can also be arranged in various lecture layouts (see dotted lines).

Worksurface

Molded fiberglass reinforced polyester top is standard. Epoxy resin top is optional. Both surfaces include integral raised perimeter edge to contain spillage.

Sink Basin

One-piece molded fiberglass polyester sink is standard. Drop-in style (epoxy resin) is optional.

Services

Unicast fixtures (cold water & double gas) are standard. Unimix fixtures (cold water, hot water and gas) are optional. Internal piping, wiring, conduit and drain lines are not included.

Connectivity

Three (3) GFI Duplex receptacles

Support Understructure

Molded one-piece fiberglass reinforced polyester understructure.

Steel Pedestal Base

Height adjustable at time of installation. Includes removable steel panels for access to services. Standard color is Black.


Accessibility

Accessibility worksurface available at ADA height (max. 34" high). Includes unicast fixtures with ADA wristblade handles.

Color Options-Epoxy

Available in Black Onyx, Lunar White, Gray, Tan, and Dark Khaki.

Color Options-Fiberglass*

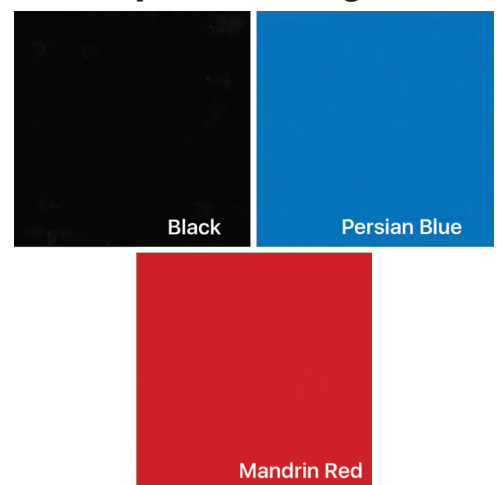

***If color match is critical, samples are available upon request.**

sales@sheldonlabs.com



800.531.7604